

Chemical Resistances for Beckman Coulter Centrifugation Products

**CAUTION**

The information provided here is from references, from current literature, or from research done by Beckman Coulter, Inc., and is only a guide for the selection of materials. No guarantee of safety based on these recommendations is expressed or implied. Many of the chemicals are explosive when concentrated or dry, or are toxic, caustic, allergenic, or carcinogenic. Observe proper handling as outlined by your laboratory safety officer.

INTRODUCTION

This table indicates the general chemical resistances of various materials to a number of chemicals commonly used in procedures involving Beckman Coulter centrifuges and accessories. You can select a chemical and determine the resistance of listed materials to the selected chemical—either satisfactory (S), marginally satisfactory (M), unsatisfactory (U), or unknown (X). Marginal resistance listings may be a combination of S and U resistances.

Materials are listed alphabetically. Chemicals are listed alphabetically by their most common name within seven categories (acids, bases, salts, gradient-forming materials, solvents, detergents, and other). Where applicable, an IUPAC (International Union of Pure and Applied Chemistry) name is shown beneath a trivial chemical name. Chemicals are either undiluted liquids or saturated (unless otherwise noted) aqueous solutions. Materials that have unsatisfactory or marginal resistance to the high concentrations used for these tests may be usable in very low (that is, millimolar) concentrations.

TEST YOUR SOLUTION UNDER OPERATING CONDITIONS IF MATERIAL PERFORMANCE IS UNCERTAIN.

Soak tests at $1 \times g$ (at 20°C) established the data for most of the materials. In some cases the resistances of tube materials also reflect their performance under centrifugation. Thus, although alcohols (for example) may be stored satisfactorily in polycarbonate or in Ultra-Clear containers, ethanol will destroy these tubes in a short period of high-speed centrifugation. This combination of material and chemical is therefore listed as U in the table. Not all combinations have been tested under the stress of centrifugation, however. Again, pretesting under actual run conditions is strongly advised.

The data for centrifuge and rotor finishes is derived mainly from splash tests in which the finish was exposed to the chemical for a matter of minutes. Satisfactory resistance under long-term exposure should not be assumed.

DECONTAMINATION OF ALUMINUM ROTORS AND ACCESSORIES

While a number of solutions are commercially marketed for use in removing radioactivity from contaminated materials, many are too harsh for use on anodized aluminum. Beckman Coulter has tested a number of solutions and found two that do not harm anodized aluminum:

- IsoClean Solution (for soaking) or RadCon Surface Spray (In U.S.A., contact Nuclear Associates [New York]; in Eastern Europe and Commonwealth States, contact Victoreen GmbH [Munich]; in South Pacific, contact Gammasonics Pty. Ltd. [Australia]; in Japan, contact Toyo Medic Co. Ltd. [Tokyo].)
- Radiacwash (In U.S.A., contact Biodex Medical Systems [Shirley, New York]; internationally, contact the U.S. office to find the dealer nearest you.)

While Beckman Coulter has tested these methods and found that they do not damage components, no guarantee of decontamination is expressed or implied. Consult your laboratory safety officer regarding the proper decontamination methods to use.

If a rotor and/or accessories are contaminated with toxic or pathogenic solutions, follow appropriate sterilization or disinfection procedures as outlined by your laboratory safety officer.

REGISTERED TRADEMARKS

No part of this document may be reproduced or transmitted in any form or by any means, electronic, mechanical, photocopying, recording or otherwise, without prior written permission from Beckman Coulter, Inc.

Beckman Coulter and the stylized logo are trademarks of Beckman Coulter, Inc. and are registered in the USPTO.

All other trademarks, service marks, products, or services are the trademarks or registered trademarks of their respective owners.

Chemicals
IUPAC Name

Chemicals IUPAC Name	acetal copolymer (celcon)	acrylic (plexiglass)	alumina (Al ₂ O ₃)	aluminum ⁷	anodic coating ⁷	nitrile rubber (Buna N)	acetal homopolymer (Delrin)	synthetic rubber (Dextrin)	epoxy resin	epoxy resin/carbon composite (Mylar)	neoprene	polyethylene terephthalate polyester film	nylon (6, 6/6)	paint, water-based	polyethylene oxide (PPO) (Noryl)	polyethylene terephthalate	polycarbonate	polyethylene terephthalate (PET)	polyethylene (HDPE)	polypropylene (LDPE)	polystyrene	polysulfone	polyurethane liner	polyurethane paint	polyvinyl chloride (PVC)	teflon (PTFE)	silicone elastomer (PPSU) (Radel)	silicone rubber (RTV) (Silastic)	titanium	flexible polyvinyl chloride	polyetherimide (PEI) (Ultem)	polyethylene terephthalate glycol (PETG) (Ultra-Clear)	synthetic rubber (Viton)	polyphenylene sulfide (PPS)			
ACIDS (aq)																																					
acetic acid (5%) ethanoic acid	S	S	S	S	S	S	S	S	S	S	S	S	M ¹	S	S	S	S	S	S	S	S	S	S	M	S	S	S	M	S	S	S	U	S	S			
acetic acid (60%) ethanoic acid	U	U	S	S	S	U	S	S	S	S	S	U	S	S	U	M	M	M	M	S	M	S	U	S	S	S	M	U	S	M	M	U	M	S			
acetic acid (glacial) ethanoic acid	U	U	S	S	S	U	M	S	S	S	U	S	U	M	S	U	U	U		U	U	U	S	U	M	S	S	U	U	S	M	U	U	U	S		
boric acid	S	S	S	U	S	S	U	S	S	S	S	S	M ⁸	S	S	S	S	S	S	S	S	S	S	S	S	M	S	S	S	S	S	S	S	S	S		
chromic acid (10%)	U	S	S	U	S	U	U	U	U	U	S	U	S	U	S	S	U	U	S	M	U	U	S	U	X	S	S	U	U	S	S	S	S	S	S		
citric acid 2-hydroxy-1,2,3-propane tricarboxylic acid	S	S	S	M	S	S	U	S	S	S	S	S	M	S	S	S	M	S	M	S	S	M	S	S	S	S	S	S	S	S	S	S	S	S	S	S	
hydrochloric acid (10%)	U	S	M	U	U	M	U	S	S	S	S	S	U	S	S	U	S	S	S	S	S	S	S	S	S	S	S	U	M	S	S	S	M	S	S		
hydrochloric acid (50%)	U	U	S	U	U	M	U	M	U	U	U	M	S	U	S	U	U	S	S	M	S	S	S	M	S	M	M	U	U	M	S	U	S	X			
iodoacetic acid 2-iodoethanoic acid	S	U ¹	S	S	S	M	S ¹	S	S	S	S ¹	M	S	S ¹	U	S	S	S ¹	S	S ¹	S ¹	U	S ¹	S	X	S	M	M	S	S	M	S	M ¹	M	X		
mercaptoacetic acid 2-mercaptoethanoic acid	S	U	S	U ⁴	S	U	S	X	M	M	S	M	S	U	U	S	U	U	U	S	S	U	U	M ¹	X	S	U	U	S	S	S	S	U	S	X		
nitric acid (10%)	U	S	M	U	S	U	U	M	U	U	S	U	S	U	S	S	S	S	U	S	U	S	S ²	S	S	S	M	S	S	S	S	M	S	S			
nitric acid (50%)	U	U	S	U	S	U	U	U	U	U	U	S	U	U	U	M	M	M	S	U	U	U	U	S ²	U	S	M	U	S	S	M	M	U	S	U		
oleic acid <i>cis</i> -9-octadenoic acid	S	S	S	S	X	S	U	M	S	S	X	U	S	S	X	S	M	M	S ²	M	S	S	M	S	S	X	U	U	S	U	S	S	M	S			

S = satisfactory resistance

M = marginal resistance

U = unsatisfactory resistance

X = unknown

= Flammability hazard. The Centrifuge instrument is not designed for use with materials capable of developing flammable or explosive vapors, or with materials which could react chemically with sufficient vigor to cause a HAZARD. Do not centrifuge such materials nor handle or store them near the centrifuge. Vapors may be ignited by exposure to electrical contacts. Depending on the centrifuge type, such exposure could occur either during normal centrifugation or under failure conditions.

¹ discoloration

² below 26°C only

explosion hazard due to possible material/chemical reaction under rotor failure conditions

⁴ dilute solutions satisfactory

⁵ below 21°C only

⁶ nonaqueous

⁷ most aluminum components have anodic coating finishes

⁸ avoid high temperatures at high concentrations

⁹ nickel acetate unsatisfactory

¹⁰ vegetable oils may be marginal or unsatisfactory

Chemicals IUPAC Name	acetal copolymer (celcon)	acrylic (plexiglass)	alumina (Al ₂ O ₃)	aluminum ⁷	anodic coating ⁷	nitrile rubber (Buna N)	acetal homopolymer	synthetic rubber (Delrin)	epoxy resin	epoxy resin/carbon composite (Mylar)	neoprene	polyethylene terephthalate polyester film	nylon (6, 6/6)	paint, water-based	polyethylene oxide (PEO) (Noryl)	polyethylene	polyethylene terephthalate (PET)	polyethylene (HDPE)	polypropylene (LDPE)	polystyrene	polysulfone	polyurethane	polyurethane liner	polyurethane paint	polyvinyl chloride (PVC)	teflon (PTFE) (Rulon A)	silicone elastomer (PPSU) (Radel)	silicone rubber	stainless steel	titanium	flexible polyvinyl chloride (PVC) (Tygon)	polyetherimide (PEI) (Ultem)	polyethylene terephthalate glycol (PETG) (Ultra-Clear)	synthetic rubber (Viton)	polyphenylene sulfide (PPS)		
ACIDS (continued)																																					
oxalic acid ethanedioic acid	U	S	S	U	U	M	U	S	S	S	S	S	M	S	S	U	S	M	S	S	S	S	S	S	S	S	M	U	S	S	S	U	S	S	S		
perchloric acid (70%)	U	S	S	U	U	U	M	U	✱	S	U	M	U	S	S	U	M	M	M	U	U	U	S	M	X	S	U	U	U	S	U	U	U	S	X		
phosphoric acid mixture (10%)	U	S	M	U	U	M	U	S	S	S	S	S	U	S	S	S	S	S	S	S	M	S	S	X	S	S	S	M	U	S	S	S	S	S	S		
phosphoric acid mixture (50%)	U	S	M	U	U	U	U	S	M	M	S	U	S	U	S	S	S	S	S	S	U	S	S	X	S	S	U	U	U	S	S	S	S	S	S		
picric acid 2,4,6-trinitrophenol	S	U	S	S	S	U	S ¹	M	M	M	S	M	S	U	S	S	U	S	S	S	U	S ¹	U	U	S	S ¹	U	M	S	M ¹	X	S	S	X			
saturated fatty acids	S	U	S	S	X	S	S	M	S	S	X	S	S	S	X	S	S	M	S	S	S	S	M	S	S	S	S	M	M	S	S	X	S	S	X		
stearic acid octadecanoic acid	S	U	S	U	X	S	S	M	S	S	X	M	S	S	S	X	S	S	M	S	S	S	S	M	S	S	M	M	S	S	S	S	X	S	X		
sulfosalicylic acid 3-carboxy-4-hydroxy- benzenesulfonic acid	S	S	S	U	U	S	S	S	S	S	S	S	U	S	S	S	S	S	S	S	X	S	S	S	S	S	U	S	S	S	S	S	S	S	X		
sulfuric acid (10%)	U	S	M	U	U	M	U	M	U	U	S	S	S	U	S	S	S	S	S	S	M	S	S	S	S	S	U	U	U	S	S	S	S	S	S	S	
sulfuric acid (50%)	U	U	M	U	U	U	U	U	U	U	M	M	S	U	S	M	U	S	S	S	S	M	U	S	S	S	S	U	U	U	S	S	U	S	M		
thioglycolic acid	U	U	S	U ⁴	S	U	S	X	M	S	S	M	S	U	U	S	U	U	U	U	S	U	U	U	M ¹	X	S	U	U	S	S	S	X	U	U	X	
trichloroacetic acid trichloroethanoic acid	U	U	S	U	U	U	U	M	S	S	U	U	S	U	M	U	M	S ²	U	S	U	U	U	S	U	X	S	U	U	U	U	M	X	U	U	S	
unsaturated fatty acids	S	S ²	S	S	X	S	U	M	S	S	X	U	X	S	S	X	S	M	M	S ²	M	M	S	M	S	S	X	M	U	S	M	S	S	M	X		
BASES (aq)																																					
ammonium hydroxide (10%)	S	S	S	U	U	S	U	S	S	S	U	S	S	U	S	U	S	S	S	S	S	S	S	S	S	S	S	S	S	S	M	M	U	S	X		
Imidazole, Amine	X	X	X	X	X	X	X	X	X	X	X	X	X	X	U	S	S	S	X	X	X	X	X	X	X	X	X	X	X	X	X	X	X	X	X	X	M

S = satisfactory resistance

M = marginal resistance

U = unsatisfactory resistance

X = unknown

☹ = Flammability hazard. The Centrifuge instrument is not designed for use with materials capable of developing flammable or explosive vapors, or with materials which could react chemically with sufficient vigor to cause a HAZARD. Do not centrifuge such materials nor handle or store them near the centrifuge. Vapors may be ignited by exposure to electrical contacts. Depending on the centrifuge type, such exposure could occur either during normal centrifugation or under failure conditions.

¹ discoloration

² below 26°C only

✱ explosion hazard due to possible material/chemical reaction under rotor failure conditions

⁴ dilute solutions satisfactory

⁵ below 21°C only

⁶ nonaqueous

⁷ most aluminum components have anodic coating finishes

⁸ avoid high temperatures at high concentrations

⁹ nickel acetate unsatisfactory

¹⁰ vegetable oils may be marginal or unsatisfactory

Chemicals
IUPAC Name

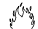
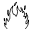


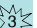
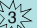
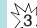
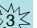
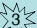

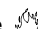
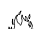
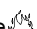
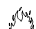
Chemicals IUPAC Name	acetal copolymer (celcon)	acrylic (plexiglass)	alumina (Al ₂ O ₃)	aluminum ⁷	anodic coating ⁷	nitrile rubber (Buna N)	acetal homopolymer (Delrin)	synthetic rubber (Dextrin)	epoxy resin	epoxy resin/carbon composite (Mylar)	polyethylene terephthalate polyester film	neoprene	polyphenylene oxide (PPO) (Noryl)	nylon (6, 6/6)	paint, water-based	polyethylene terephthalate (PET)	polycarbonate	polyethylene terephthalate (PET)	polyethylene (HDPE)	polypropylene	polystyrene	Polysulfone	polyurethane liner	polyurethane paint	polyvinyl chloride (PVC)	teflon (PTFE)	polyphenylsulfone (PPSU) (Radcl)	silicone elastomer (RTV) (Silastic)	stainless steel	titanium	flexible polyvinyl chloride	polyetherimide (PEI) (Ultem)	polyethylene terephthalate glycol (PETG) (Ultra-Clear)	synthetic rubber (Viton)	polyphenylene sulfide (PPS)		
SALTS (aq) (continued)	S	S	S	M	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S
magnesium chloride	S	S	S	M	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	
nickel salts	S	S	S	U	S	S	S	S	S	S	M	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	M	S	S	S	S	S	S ⁹	S	S	
potassium bromide	S	S	S	U	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	M	S	S	S	S	S	S	S	S	S	
potassium carbonate	S ²	M ¹	S	M	U	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	
potassium chloride	S	S	M	U	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	U	S	S	S	S	S	S	S	S	S	
potassium permanganate	S	S ²	S ¹	S	S	S	S	S	S	S	S ⁴	S	U	S	S	S	S	S	S ²	M	S ¹	U	M	U	S	S	S ¹	S ¹	M	S	U	S	S ¹	M	S	S	
silver nitrate	S	S	S	U	S	M	S	S	S	S	S ⁴	S	S	S	S	S	S	S	S	S	S	X	M	S	S	S	S	M	S	S	S	S	S	S	S	S	S
sodium borate	S	S	S	S	S	S	S	S	S	S	S	S	U	S	S	S	S	S	S	S	S	S	S	S	S	M	S	M	S	S	S	S	S	S	S	S	
sodium carbonate	S	S	S	M	U	S	S	S	S	S	S	S	S	S	S	S	S	S	M	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	
sodium chloride	S	S	S	U	U	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	M	S	S	S	S	S	S	S	S	S	S
sodium nitrate	S	S	S	S	S	S	S	S	S	S	M	S	S	S	U	S	S	S	S	S	M	S	S	S	S	U	S	S	S	S	S	S	S	S	S	S	S
sodium sulfate	S	S	S	U	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	M	S	S	S	S	S	M	S	S	S	S	S	S	S	S	S	S
sodium sulfite	S	S	S	S ²	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	M	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	X
zinc chloride	S	S	S	U	U	S	U	S	S	S	S	S	S	S	S	S	S	S	S	S	U	S	S	S	S	S	U	S	S	S	S	S	S	S	S	S	S
zinc sulfate	S	S	S	U	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	U	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S
GRADIENT FORMING MATERIALS (aq)																																					
cesium acetate	S	S	S	M	X	S	S	S	S	X	S	S	S	S	S	S	S	S	S	S	U	S	S	S	S	S	M	S	S	S	S	X	S	S	X		
cesium ethanoate	S	S	S	M	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	M	S	S	S	S	S	S	S	S	S	X	
cesium bromide	S	S	S	M	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	M	S	S	S	S	S	S	S	S	S	X	

S = satisfactory resistance
M = marginal resistance
U = unsatisfactory resistance
X = unknown

= Flammability hazard. The Centrifuge instrument is not designed for use with materials capable of developing flammable or explosive vapors, or with materials which could react chemically with sufficient vigor to cause a HAZARD. Do not centrifuge such materials nor handle or store them near the centrifuge. Vapors may be ignited by exposure to electrical contacts. Depending on the centrifuge type, such exposure could occur either during normal centrifugation or under failure conditions.

¹ discoloration
² below 26°C only
³ explosion hazard due to possible material/chemical reaction under rotor failure conditions
⁴ dilute solutions satisfactory
⁵ below 21°C only
⁶ nonaqueous
⁷ most aluminum components have anodic coating finishes
⁸ avoid high temperatures at high concentrations
⁹ nickel acetate unsatisfactory
¹⁰ vegetable oils may be marginal or unsatisfactory

Chemicals
IUPAC Name

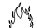
	acetal copolymer (celcon)	acrylic (plexiglass)	alumina (Al ₂ O ₃)	aluminum ⁷	anodic coating ⁷	nitrile rubber (Buna N)	acetal homopolymer (Delrin)	synthetic rubber (EPDM)	epoxy resin	epoxy resin/carbon composite (Mylar)	polyethylene terephthalate polyester film	neoprene	polyphenylene oxide (PPO) (Noryl)	nylon (6, 6/6)	paint, water-based	polyethylene terephthalate (PET)	polycarbonate	polyethylene (HDPE)	polyethylene (LDPE)	polypropylene	polystyrene	polysulfone	polyurethane liner	polyurethane paint	polyvinyl chloride (PVC)	polyphenylsulfone (PPSU) (Radel)	teflon (PTFE) (Rulon A)	silicone elastomer (RTV) (Silastic)	stainless steel	titanium	flexible polyvinyl chloride (PVC) (Tygon)	polyetherimide (PEI) (Ultem)	polyethylene terephthalate glycol (PETG) (Ultra-Clear)	synthetic rubber (Viton)	polyphenylene sulfide (PPS)		
SOLVENTS																																					
acetone 2-propanone 	M	U	S	M	S	U	S	S	U	U	M	U	U	S	U	M	U	S	S	S	U	U	U	M	U	M	S	U	M	M	S	U	S	U	U	S	
acetonitrile ethanenitrile 	S	U	S	S	S	U	S	U	M	M	S	S	U	S	U	S	S	S	S	U	U	U	U	M	M	S	S	S	S	S	U	X	U	U	X		
benzene 	M	U	S	S	S	U	M	U	U	U	S	U	U	S	U	U	U	U	U	U	U	U	S	U	M	S	U	U	M	S	U	X	U	S	S		
carbon tetrachloride tetrachloromethane	S	U	S			M	S	U	U	U	S	U	U	S	S	S	U	U	U	U	U	S	U	S	M	U	U	U		U	S	U	S	S			
chloroform trichloromethane	S	U	S			U	S	U	S	S	S	U	U	U	S	S	U	U	U	U	U	U	M	U	U	S	U	U	U		U	U	U	U	M		
cresol mixture methylphenol	S	U	S	S	S	U	S	U	U	U	S	U	U	U	S	U	U	U	U	U	U	U	U	X	S	S	U	S	S	U	X	U	S	S			
cyclohexane 	S	U	S	S	S	S	U	S	S	S	U	U	S	S	S	U	U	U	U	U	M	U	S	S	U	U	M	M	U	S	U	S	U	S			
diethyl ether ethoxyethane 	S	U	S	S	S	U	S	U	S	S ²	S	U	U	S	S	S	U	U	U	U	U	U	S	U	S	S	U	S	U	S	U	S	U	S	U	S	
diethyl ketone 3-pentanone 	S	U	S	S	X	U	S	M	M	M	X	U	U	S	U	X	U	U	U	M	U	U	S	U	M	S	X	U	M	S	U	S	U	U	X		
N,N-dimethylformamide N,N-dimethylmethanamide 	S	U	S	S	S	M	S	M	M	M	S	U	U	S	S	S	U	S	S	S	U	U	U	S	U	X	S	S	M	S	S	U	X	U	U	S	
dimethyl sulfoxide sulfanylbis[methane]	S	M	S	S	S	U	S	S	S	S	U	S	S	S	S	U	S	S	S	S	U	U	S	U	X	S	U	S	S	S	U	X	U	U	X		
dioxane 1,4-dioxacyclohexane 	S	U	S	M	S	U	M	M	S	S	S	U	U	S	U	S	U	U	U	S	U	U	U	S	U	X	S	S	U	S	S	U	X	U	U	X	

S = satisfactory resistance

M = marginal resistance


U = unsatisfactory resistance

X = unknown

 = Flammability hazard. The Centrifuge instrument is not designed for use with materials capable of developing flammable or explosive vapors, or with materials which could react chemically with sufficient vigor to cause a HAZARD. Do not centrifuge such materials nor handle or store them near the centrifuge. Vapors may be ignited by exposure to electrical contacts. Depending on the centrifuge type, such exposure could occur either during normal centrifugation or under failure conditions.

¹ discoloration

² below 26°C only

 explosion hazard due to possible material/chemical reaction under rotor failure conditions

⁴ dilute solutions satisfactory

⁵ below 21°C only

⁶ nonaqueous

⁷ most aluminum components have anodic coating finishes

⁸ avoid high temperatures at high concentrations

⁹ nickel acetate unsatisfactory

¹⁰ vegetable oils may be marginal or unsatisfactory

Chemicals
IUPAC Name

acetal copolymer (celcon)
acrylic (plexiglass)
alumina (Al₂O₃)
aluminum⁷
anodic coating⁷
nitrile rubber (Buna N)
acetal homopolymer (Delrin)
synthetic rubber (EPRM)
epoxy resin
epoxy resin/carbon composite (Mylar)
polyethylene terephthalate polyester film
neoprene
polyphenylene oxide (PPO) (Noryl)
nylon (6, 6/6)
paint, water-based
polyethylene-based
polycarbonate
polyethylene terephthalate (PET)
polyethylene (HDPE)
polypropylene
polystyrene
polysulfone
polyurethane liner
polyurethane paint
polyvinyl chloride (PVC)
teflon (PTFE)
polyphenylsulfone (PPSU) (Radel)
silicone elastomer (RTV) (Silastic)
silicone rubber
stainless steel
titanium
flexible polyvinyl chloride (PVC) (Tygon)
polyetherimide (PEI) (Utem)
polyethylene terephthalate glycol (PETG) (Ultra-Clear)
synthetic rubber (Viton)
polyphenylene sulfide (PPS)

SOLVENTS (continued)	S	U	S	S	S	S	S	S	S	S	S	S	S	S	S	U	S	S	S	M	S	S	S	S	S	M	S	M	S	M	S	M	S	U	S	S	
ethanol (50%)	S	U	S	S	S	S	S	S	S	S	S	S	S	S	S	U	S ⁵	S	S	M	S	S	S	S	S	M	S	M	S	M	S	M	S	U	S	S	
ethanol (95%)	S	U	S	S	S	S	S	S	S	S	S	S	S	S	S	U	S ⁵	S	S	M	S	S	S	S	S	M	S	M	S	M	S	M	S	U	S	S	
ether	S	U	S	S	S	U	S	U	S	S ²	S	U	U	S	S	S	U	U	U	U	S	U	U	S	U	U	S	M	S	M	S	M	S	U	U	S	
ethyl acetate ethyl ethanoate	S	U	S	M	S	U	S	M	M	M	S ¹	U	U	S	U	U	U	S	U	U	U	S	U	U	S	U	M	M	S	U	S	U	U	S			
ethylene glycol 1,2-ethanediol	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	U	S	M	S	S	S	S	S	S	S	S	M	S	M	S	M ⁵	S	S	S	S		
hexane	S	S	S	S	S	S	S	U	S	S ²	S	S	U	S	S	S	U	U	U	M ¹	U	S	S	S	U	S	S	U	U	S	S	U	S	U	S	S	
isopropyl alcohol 2-propanol	S	U	S	M	S	S	S	S	S	S ²	S	M	S	S	S	S	U	S	S	S	M	S	S	S	S	S	S	S	M	S	M	S	M	S	S	S	
kerosene	S	S	S	S	S	S	S	U	S	S	S	U		S	S	S	U	U	U	U	U	S	S	S	S	S	S	U	U	S	S	U	S	U	S	S	
methanol	S	U	S	S	S	S	S	S	S	S	S	S		S	S	S	S	U	S	M	S	M	M	S	S	S	M	S	M	S	M	S	M	S	U	M	S
methylene chloride dichloromethane	M	U	S			U	S	U	U	S	S	U	U	U	U	S	U	U	U	U	U	U	U	S	U	S	U	U	M		U	U	U	M	M		
methyl ethyl ketone 2-butanone	S	U	S	S	S	U	S	S	M	M	S	U	U	S	U	S	U	U	S ²	U	U	U	M	U	M	S	U	U	S	U	S	U	S	U	U	M	
phenol (5%)	U	U	S	U	S	U	M	M	M	M	S	U	M	U	M	S	U	U	U	S	M	U	U	M	U	U	S	U	U	U	S	M	S	U	S	S	
phenol (50%)	U	U	S	U	S	U	U	U	U	U	U	M	U	U	U	U	U	U	S ²	U	U	U	M	U	U	S	U	U	U	S	M	S	U	S	S ¹¹		
tetrahydrofuran	M	U	S	S	S	U	U	U	U	U	S	U	U	S	U	S	U	U	U	U	U	U	U	U	U	S	S	U	U	S	S	U	X	U	U	S	
toluene methylbenzene	S	U	S	S	S	U	S	U	M	M	S	U	U	S	S	S	U	U	U	U	U	U	U	U	U	M	U	M	S	U	U	S	S	U	S	S	

S = satisfactory resistance

M = marginal resistance

U = unsatisfactory resistance

X = unknown

= Flammability hazard. The Centrifuge instrument is not designed for use with materials capable of developing flammable or explosive vapors, or with materials which could react chemically with sufficient vigor to cause a HAZARD. Do not centrifuge such materials nor handle or store them near the centrifuge. Vapors may be ignited by exposure to electrical contacts. Depending on the centrifuge type, such exposure could occur either during normal centrifugation or under failure conditions.

¹ discoloration

² below 26°C only

= explosion hazard due to possible material/chemical reaction under rotor failure conditions

= dilute solutions satisfactory

⁵ below 21°C only

⁶ nonaqueous

⁷ most aluminum components have anodic coating finishes

⁸ avoid high temperatures at high concentrations




⁹ nickel acetate unsatisfactory

¹⁰ vegetable oils may be marginal or unsatisfactory

¹¹ satisfactory resistance when exposed to phenol (70%)

Chemicals
IUPAC Name

acetal copolymer (Celcon)
acrylic (plexiglass)
alumina (Al₂O₃)
aluminum⁷
anodic coating⁷
nitrile rubber (Buna N)
acetal homopolymer (Dextrin)
synthetic rubber (EPDM)
epoxy resin
epoxy resin/carbon composite (Mylar)
polyethylene terephthalate
neoprene
polyphenylene oxide (PPO) (Noryl)
Paint, water-based
polyethylene terephthalate (PET)
polycarbonate
polyethylene (HDPE)
polyethylene (LDPE)
polystyrene
polysulfone
polyurethane liner
polyurethane paint
polyvinyl chloride (PVC)
teflon (PTFE) (Rulon A)
silicone elastomer (PPSU) (Radel)
silicone rubber
stainless steel
titanium
flexible polyvinyl chloride (PVC) (Tygon)
polyetherimide (PEI) (Ultem)
polyethylene terephthalate glycol (PETG) (Ultra-Clear)
synthetic rubber (Viton)
polyphenylene sulfide (PPS)


Chemicals IUPAC Name	acetal copolymer (Celcon)	acrylic (plexiglass)	alumina (Al ₂ O ₃)	aluminum ⁷	anodic coating ⁷	nitrile rubber (Buna N)	acetal homopolymer (Dextrin)	synthetic rubber (EPDM)	epoxy resin	epoxy resin/carbon composite (Mylar)	polyethylene terephthalate	neoprene	polyphenylene oxide (PPO) (Noryl)	Paint, water-based	polyethylene terephthalate (PET)	polycarbonate	polyethylene (HDPE)	polyethylene (LDPE)	polystyrene	polysulfone	polyurethane liner	polyurethane paint	polyvinyl chloride (PVC)	teflon (PTFE) (Rulon A)	silicone elastomer (PPSU) (Radel)	silicone rubber	stainless steel	titanium	flexible polyvinyl chloride (PVC) (Tygon)	polyetherimide (PEI) (Ultem)	polyethylene terephthalate glycol (PETG) (Ultra-Clear)	synthetic rubber (Viton)	polyphenylene sulfide (PPS)			
DETERGENTS (continued)																																				
Vesphen 11	S	S	X	X	X	S	S	S	S	S	S	S	S	X	S	U	S	M	S	U	S	X	S	S	S	S	X	X	S	S	M	S	S			
Wescodyne	S	S	X	X	X	S	S	S	S	S	S	S	S	X	S	U	S	M	S	U	S	X	S	S	S	S	X	X	S	S	M	S	S			
Zephiran chloride (1%)	S	S	S	M	S	S	S	S	S	S	S	S	S	S	S	S	M	S	S	S	X	S	S	S	S	S	S	S	U	S	M	S	S			
OTHER																																				
n-butyl phthalate⁶ dibutyl 1,2-benzenedicarboxylate	S	U	S	S	S	U	S	S	X	S	S	U	U	S	S	S	U	U	U	U	U	S	S	U	M	M	S	U	X	M	S	M				
dibutyl phthalate	S	S	S	S	S	U	S	S	X	S	S	U	U	S	S	U	U	U	U	U	S	X	S	U	S	M	M	S	U	S	X	S	X			
diethyl pyrocarbonate pyrocarbonic acid diethyl ester	S	M	S	S	S	U	S	X	S	S	S	U	S	S	U	S	S	S	U	S	X	U	M ¹	S	S	S	S	S	S	S	U	S	X			
ethylene oxide vapor⁶ oxirane 	X	S	X	S	X	U	X	U	U	U	X	U	X	S	M	S	S	S	S	S	U	S	X	S	U	S	S	S	U	U	S	S	S	U	U	U
formaldehyde methanal	S	S	S	S	S	S	S	S	S	S	S	S	S	X	S	S	S	S	S	U	S	U	S	S	S	S	M	S	M	S	S	S	S			
formalin (40%)	X	S	S	S	S	S	X	X	S	S	S	U	S	S	S	S	S	S	M	S	U	S	U	S	S	S	M	M	S	M	S	S	S	X		
hydrogen peroxide (3%)	S	S	S	S	U	S	S	S	U	U	S	S ²	S	U	S	S	S	S	S	S	S	M	S	S	S	S	S	S	S	S	S	S	S	S	S	
hydrogen peroxide (10%) 	U	S	S	U	U	S	U	U	U	U	S	S ²	S	U	S	S	S	S	S	S	S	U	S	S	S	S	M	S	U	S	S	S	S	S		
2-mercaptoethanol 	S	U	S	S	S	U	S	S	M	M	S	U	U	S	U	S	S	S	S	S	S	X	S	U	S	S	S	S	S	S	S	U	S	X		

S = satisfactory resistance

M = marginal resistance

U = unsatisfactory resistance

X = unknown

 = Flammability hazard. The Centrifuge instrument is not designed for use with materials capable of developing flammable or explosive vapors, or with materials which could react chemically with sufficient vigor to cause a HAZARD. Do not centrifuge such materials nor handle or store them near the centrifuge. Vapors may be ignited by exposure to electrical contacts. Depending on the centrifuge type, such exposure could occur either during normal centrifugation or under failure conditions.

¹ discoloration

² below 26°C only

 explosion hazard due to possible material/chemical reaction under rotor failure conditions

⁴ dilute solutions satisfactory

⁵ below 21°C only

⁶ nonaqueous

⁷ most aluminum components have anodic coating finishes

⁸ avoid high temperatures at high concentrations

⁹ nickel acetate unsatisfactory

¹⁰ vegetable oils may be marginal or unsatisfactory

Chemicals
IUPAC Name

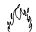
Chemicals IUPAC Name	acetal copolymer (celcon)	acrylic (plexiglass)	alumina (Al ₂ O ₃)	aluminum ⁷	anodic coating ⁷	nitrile rubber (Buna N)	acetal homopolymer	synthetic rubber (Delrin)	epoxy resin	epoxy resin/carbon composite (Mylar)	neoprene	polyethylene terephthalate polyester film	nylon (6, 6/6)	paint, water-based	polyethylene oxide (PPO) (Noryl)	polycarbonate	polyethylene terephthalate (PET)	polyethylene (HDPE)	polyethylene (LDPE)	polystyrene	polysulfone	polyurethane liner	polyurethane paint	polyvinyl chloride (PVC)	teflon (PTFE)	polyphenylsulfone (PPSU) (Radcl)	silicone elastomer (RTV) (Radcl)	stainless steel	titanium	flexible polyvinyl chloride (PVC) (Tygon)	polyetherimide (PEI) (Ultem)	polyethylene terephthalate glycol (PETG) (Ultra-Clear)	synthetic rubber (Viton)	polyphenylene sulfide (PPS)		
OTHER (continued)	S	S	S	S	X	S	S	U	S	S	S	M	S	S	S	S	M	U	S	U	U	S	S	S	S	M	U	S	S	M	S	S	S	S		
oils (petroleum)	S	S	S	S	X	S	S	U	S	S	S	M	S	S	S	M	U	S	U	U	S	U	S	S	S	X	S	S	S	S	M	S	S	S		
oils (other)	S	S ¹⁰	S	S	S	S	S	U	S	S	S	U	S ¹⁰	S	S	S	M	U	U	S	U	S	S	X	S	S	S	M	S	S	M	S	M	S	S	
physiologic media (e.g., culture media, milk, serum, urine)	S	S	S	M	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	M	S	S	S	S	S	S	S	X	
sodium hypochlorite (5 1/4% solution; unscented commercial bleach)	U	S	M	U	S	M	U	S	S	S	S	U	S	U	S	S	S	S	S	M ⁴	S	S	U	M	S	S	S	S	M	U	S	M	S	S	S	
Tris buffer (neutral pH) tris (hydroxymethyl) aminomethane	S	S	S	U	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	S	X
urea	S	S	S	S	S	U	S	S	X	X	S	S	S	S	S	S	S	S	S	S	M	X	S	S	S	S	S	M	S	S	S	S	S	S	S	S

S = satisfactory resistance

M = marginal resistance


U = unsatisfactory resistance

X = unknown

 = Flammability hazard. The Centrifuge instrument is not designed for use with materials capable of developing flammable or explosive vapors, or with materials which could react chemically with sufficient vigor to cause a HAZARD. Do not centrifuge such materials nor handle or store them near the centrifuge. Vapors may be ignited by exposure to electrical contacts. Depending on the centrifuge type, such exposure could occur either during normal centrifugation or under failure conditions.

¹ discoloration

² below 26°C only

 explosion hazard due to possible material/chemical reaction under rotor failure conditions

⁴ dilute solutions satisfactory

⁵ below 21°C only

⁶ nonaqueous

⁷ most aluminum components have anodic coating finishes

⁸ avoid high temperatures at high concentrations

⁹ nickel acetate unsatisfactory

¹⁰ vegetable oils may be marginal or unsatisfactory

GLOSSARY OF TERMS

anodized coating	a thin, hard layer of aluminum oxide formed electrochemically on aluminum rotor and/or accessory surfaces as a protective coating for corrosion resistance
Buna N	black nitrile rubber used for O-rings and gaskets in rotor assemblies
Delrin	thermoplastic material (acetal homopolymer) used for most tube adapters
EPDM	ethylene propylene rubber used for O-rings and pad adapters
HDPE	high density polyethylene used for adapters
LDPE	low density polyethylene used for tubes and bottles
neoprene	black synthetic elastomer used for O-rings in some tube caps and bottle cap assemblies
Noryl	modified thermoplastic polyphenylene oxide (PPO) used for floating spacers (part of the <i>g</i> -Max system) and some polycarbonate bottle caps
PET	polyethylene terephthalate used in some adapters
Polyphenylene Sulfide	(PPS) used in cap closures and other products as specified
Radel	polyphenylsulfone (PPSU) used in plugs, cap closures, cannisters, and other accessories
Ultem	polyetherimide (PEI) used in adapters, covers, and spacers
Viton	fluorocarbon elastomer used in high-temperature applications

Trademarks

Avanti, Beckman Coulter, and the stylized logo are trademarks of Beckman Coulter, Inc. and are registered in the USPTO. All other trademarks, service marks, products, or services are trademarks of their respective holders.



Beckman Coulter, Inc. • 250 S. Kraemer Blvd. • Brea, California 92821 U.S.A.
Sales and Service: 1-800-742-2345 • Internet: www.beckmancoulter.com • Printed in U.S.A.

© 2014 Beckman Coulter, Inc.
All rights reserved